



DEVICE CONNECTIVITY FOR
POINT-OF-CARE TESTING

RALS™ DEVICE INTERFACING

For rapid whole blood testing from the Cardiovascular Operating Room to the Cardiac Catheterization Lab, Hemochron Signature Elite system standardizes ACT+, ACT-LR, PT and APTT testing at the point-of-care. Rapid, accurate and reproducible results, for patient management decisions in acute care settings.

Connectivity by the RALS System enables users to:


- Centralize all POC test results into the RALS System
- Transfer results to an EMR/EHR/LIS
- Email results to operators and managers
- Access and manage mis-identified patient and out of range QC results
- Monitor and manage operator certifications



VISIT RALS.COM

RALS™ | Hemochron® Signature Elite Connectivity

HOSPITAL / FACILITY / INSTITUTION NETWORK

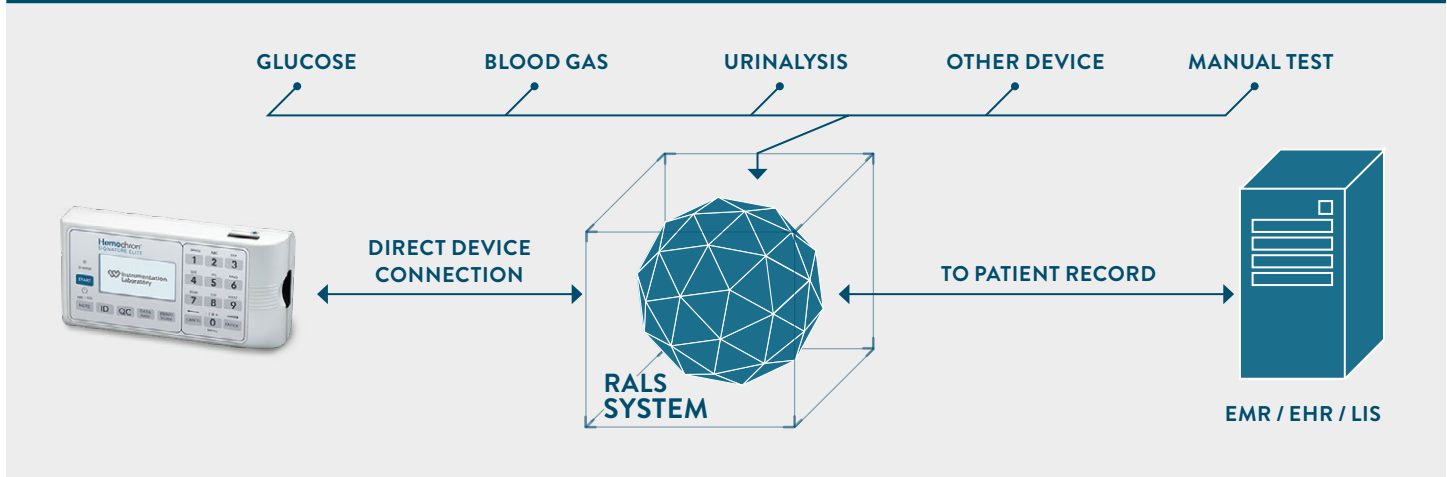

DATA FLOW

The Hemochron Signature Elite transfers data to RALS via a dataport or RALS Remote Connect software, where it is evaluated and either sent directly to the EMR/EHR/LIS or held for review.

BI-DIRECTIONAL COMMUNICATION

From device to RALS:

- Patient and QC Results
- Patient ID and Operator ID
- Date/Time
- Location

From RALS to device:

- Operator ID
- Reagents/Tests Materials
- Device Configuration

AVAILABLE REPORTS FROM THE RALS SYSTEM

- Anomalous Results by Device
- Anomalous Results by Operator
- Institution Statistics
- Location Statistics
- LIS Statistics by Location
- LIS Statistics by Operator
- Operator ID Compliance by Location
- Operator Pending Recertification Date
- Operator Pending Recertification Range
- Operator Statistics
- Patient ID Compliance by Operator
- Patient ID Compliance by Location
- Patient Results by Operator
- Patient Results by Patient
- QC Test Results
- Device Statistics Summary
- Levey-Jennings
- Corrective Action Compliance by Location and by Operator
- Patient Result Comments
- Comment Codes
- Test Materials
- Workload Summary
- Operators with No Certifications
- RALS Device Log

HOW THE DEVICE CONNECTS

- RALS Remote Connect
- RALS Dataport

TO LEARN MORE, VISIT RALS.COM

© 2019 Abbott. All rights reserved. All trademarks referenced are trademarks of either the Abbott group of companies or their respective owners. Any photos displayed are for illustrative purposes only. Any person depicted in such photos is a model. INFO216904 AI-COLL-0052 Rev 2 06/19

