



DEVICE CONNECTIVITY  
FOR POINT-OF-CARE TESTING

# RALS™ DEVICE INTERFACING

The Medtronic ACT Plus® System is a microprocessor controlled instrument designed to determine coagulation endpoints in fresh and citrated whole blood samples.


*Connectivity by the RALS System enables users to:*

- Centralize all POC test results into the RALS System
- Transfer results to an EMR/EHR/LIS
- Email results to operators and managers
- Access mis-identified patient and out of range results
- Monitor and manage operator certifications



VISIT [RALS.COM](http://RALS.COM)

# RALS™ | Medtronic ACT Plus® Connectivity

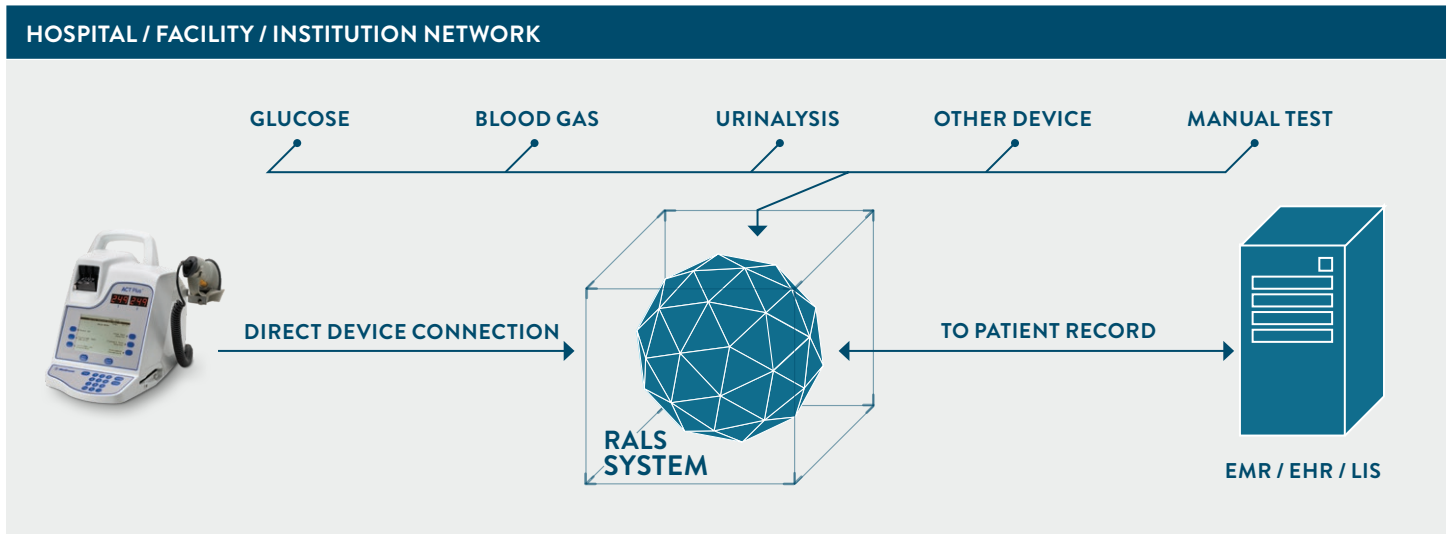

## DATA FLOW

The Medtronic ACT Plus System transfers data to RALS via RALS Remote Connect software or a RALS Dataport, where it is evaluated and either sent directly to the EMR/EHR/LIS or held for review.

## DIRECT-TO-DEVICE INTERFACING

The Medtronic ACT Plus System transfers a wide range of data to the RALS System, including:

- Patient test results
- QC test results
- EQC results
- Operator ID
- Date/Time
- Location

## UNI-DIRECTIONAL COMMUNICATION

From device to the RALS System:

- Patient and QC Results
- Patient ID and Operator ID
- Date/Time
- Location

## AVAILABLE REPORTS

- Anomalous Results by Device
- Anomalous Results by Operator
- Institution Statistics
- Location Statistics
- LIS Statistics by Location
- LIS Statistics by Operator
- Operator ID Compliance by Location
- Operator Pending Recertification Date
- Operator Pending Recertification Range
- Operator Statistics
- Patient ID Compliance by Operator
- Patient ID Compliance by Location
- Patient Results by Operator
- Patient Results by Patient
- QC Test Results
- Device Statistics Summary
- Corrective Action Compliance by Location and by Operator
- Comment Codes
- Test Materials
- Workload Summary
- Operators with No Certifications
- RALS Device Log

## HOW THE DEVICE CONNECTS

- RALS Remote Connect (RRC)
- RALS Dataport

TO LEARN MORE, VISIT [RALS.COM](http://RALS.COM)

© 2019 Abbott. All rights reserved. All trademarks referenced are trademarks of either the Abbott group of companies or their respective owners. Any photos displayed are for illustrative purposes only. Any person depicted in such photos is a model. AI-COLL-0054 Rev 2 05/19

