



DEVICE CONNECTIVITY FOR  
POINT-OF-CARE TESTING

# RALS™ DEVICE INTERFACING

HemoCue® Hb 201 DM systems offer the lab accuracy and ease of use that has become a standard in point-of-care testing, with features like QC lockout and built-in barcode scanning for patient ID to support patient safety and privacy.

## CONNECTIVITY TO RALS ENABLES USERS TO:


- Centralize all test results into the RALS system
- Transfer results to an EMR/EHR/LIS
- Email results to operators and managers
- Access and manage mis-identified patient and out of range QC results



VISIT [RALS.COM](http://RALS.COM)

# RALS™ | HemoCue® Hb System Interfacing

## HOSPITAL / FACILITY / INSTITUTION NETWORK

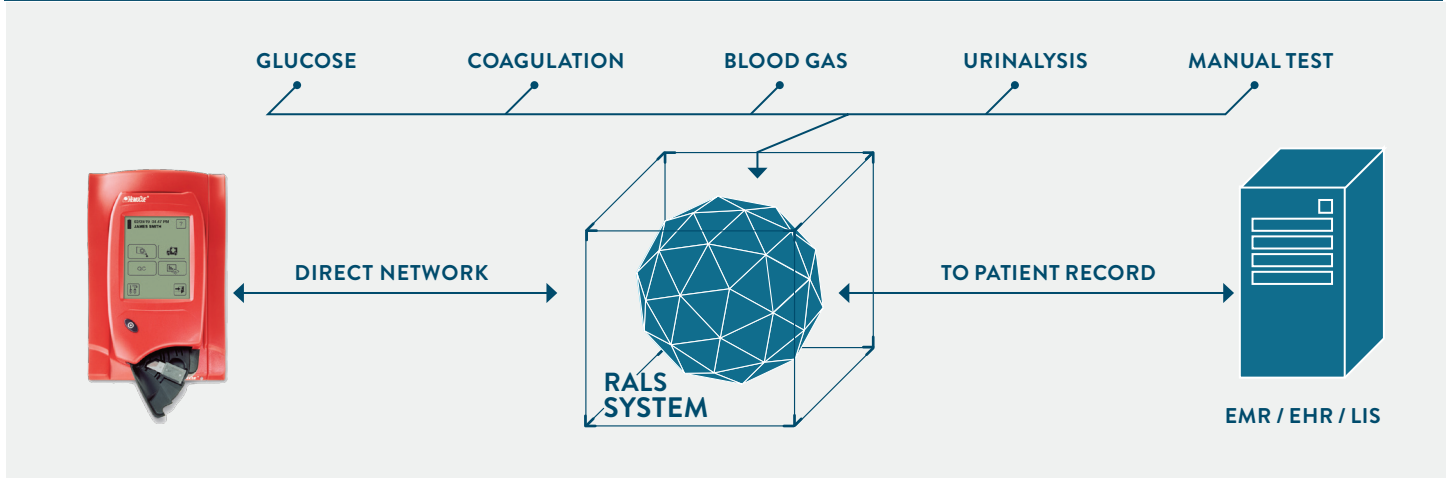

## DATA FLOW

HemoCue® Hb systems transfer data to RALS via a direct network connection, using the device's base unit, where it is evaluated and either sent directly to the EMR/EHR/LIS or held for review.

## BI-DIRECTIONAL COMMUNICATION

Data sent from the device to RALS™

- Comments (device errors)
- Materials
- Operators
- Patients
- Accession
- Results
- Log Notes as device events

Data sent from RALS™ to the device

- Operator List
- Patient List
- Test Materials
- Comments
- Instrument Configurations

## AVAILABLE REPORTS FROM THE RALS SYSTEM

- RALS™ Devices Log
- Workload Summary
- Anomalous Results by Device
- Anomalous Results by Operator
- Operator Pending Recertification Range
- Operator Pending Recertification Date
- Operators with No Certifications
- Operator Statistics
- Patient Results by Operator
- Patient Results by Patient
- LIS Statistics by Location
- LIS Statistics by Operator
- Institution Statistics
- Location Statistics
- Levey-Jennings
- QC Test Results
- Device Statistics Summary
- Test Materials

TO LEARN MORE, VISIT [RALS.COM](https://www.rals.com)

© 2019 Abbott. All rights reserved. All trademarks referenced are trademarks of either the Abbott group of companies or their respective owners. Any photos displayed are for illustrative purposes only. Any person depicted in such photos is a model. AI-COLL-0122 Rev. 1 12/19

